



240V Dog Teeth Cleaning Ultrasonic Machine USB Ultrasonic Cleaner For Dogs Teeth

Our Product Introduction

for more products please visit us on product-pet.com

Basic Information

- Place of Origin: Sichuan, China
- Brand Name: Petbelief
- Model Number: PLQURE01
- Minimum Order Quantity: 100 sets
- Price: Negotiable



Product Specification

- Feature: Sustainable
- Application: Dogs
- Type: Brushes
- Material: Plastic
- Power Source: CHARGE
- Charging Time: 4 Hour
- Voltage: 110-240V
- Product Name: Pet Plaque Remover For Teeth Kit
- MOQ: 100 Pcs
- Color: Black/White
- Battery: USB Rechargeable
- Function: Deep Clean
- Grooming Products Type: Grooming Tools
- Highlight: Dog Teeth Cleaning Ultrasonic Machine, USB Ultrasonic Cleaner For Dogs Teeth, 240V Ultrasonic dog teeth cleaner



More Images



Product Description

Teeth Pet Ultrasonic Tooth Cleaner Kit with Rechargeable Dental Care Tools

Products Description

1. Introduction

Worried about your pet's dental health? Our Teeth Pet Ultrasonic Tooth Cleaner Kit offers a comprehensive solution to keep your pet's teeth clean and healthy. Designed with advanced ultrasonic technology, this rechargeable dental care tool ensures effective plaque removal and promotes overall oral hygiene for your furry friends.

2. Key Features

- "Ultrasonic Cleaning Technology": Utilizes high-frequency vibrations to effectively remove plaque and tartar, ensuring thorough cleaning.
- "Rechargeable Battery": Equipped with a long-lasting, rechargeable battery, providing convenience and sustainability.
- "Multiple Cleaning Modes": Features different intensity settings to cater to various dental needs and sensitivities.
- "Ergonomic Design": Designed for comfortable handling, making it easy to clean your pet's teeth without causing discomfort.

3. Usage Scenarios

Perfect for regular dental care at home, this kit is ideal for pet owners who want to maintain their pet's oral hygiene without frequent vet visits. Use it to clean your pet's teeth after meals, during grooming sessions, or as part of their regular health routine. Suitable for both dogs and cats, this tool ensures your pet's teeth stay clean and healthy.

4. Safety Tips

Always start with the lowest intensity setting to ensure your pet's comfort. Gradually increase the intensity as your pet gets used to the cleaning process. Regularly check the tool for any signs of wear and tear, and replace the cleaning heads as needed. Keep the tool out of reach of children and store it in a safe place.

5. Customer Reviews

"This tooth cleaner is amazing! My dog's teeth have never looked better. The ultrasonic technology is very effective, and the rechargeable battery is a huge plus. Highly recommend for any pet owner!" - Mark T.

Produce Name	Pet Plaque Remover for Teeth
Material	Plasitc
Colors	2 Color
MOQ	100 Pcs

NOTICE This product is a patented product, with the European Union, the United States and other patents, imitators must be investigated, whoever produces and sells imitations will be held legally responsible for compensation!

Product Size



Product Details

BEFORE AND AFTER

	BEFORE	AFTER
		
		
		

How to Use



1.Remove or replace the cleaning head as you need



2.Use the wrench to tight up the cleaning head



3.Press the machine to choose the suitable mode and move the cleaning head up and down



4.Turn off the operation and use manual tools to clean the rest calculus

Ultrasonic High Efficiency Teeth Cleaning

! NOTE


This product uses ultrasonic technology and does not vibrate and make sound when turned on, but begins to vibrate when the product comes into contact with hard surfaces (e.g., teeth).



Before use



After use

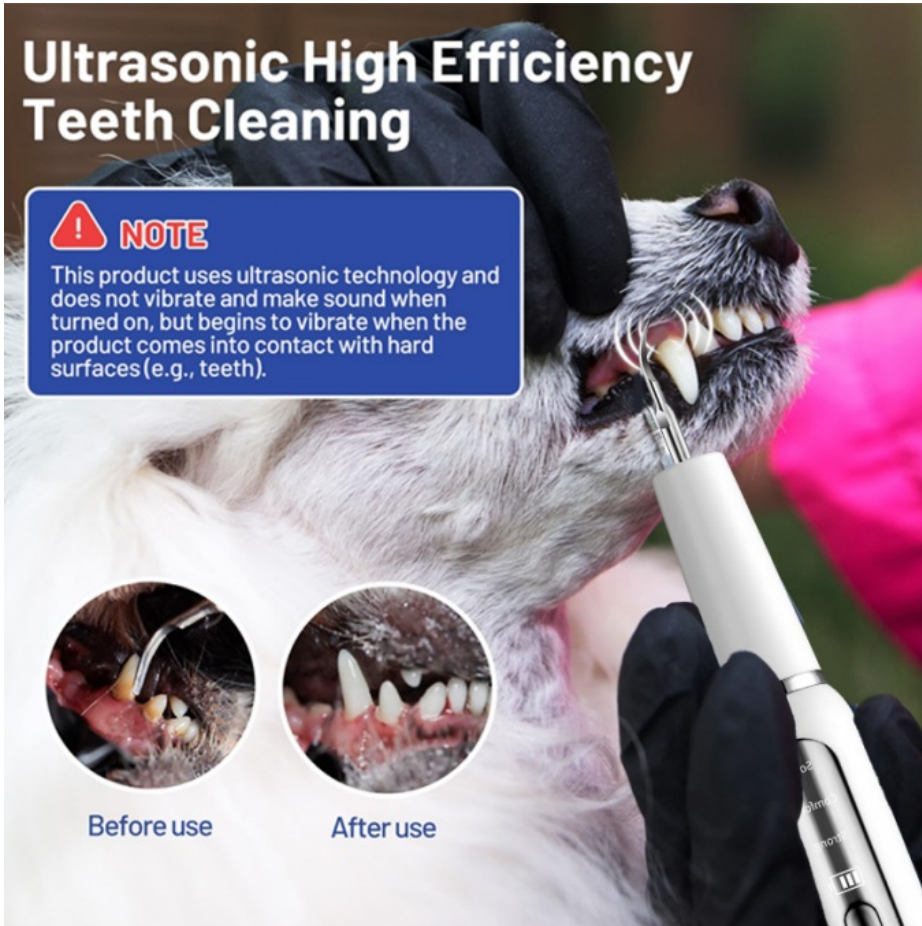

So safe for gums

Turn on NO Response! Only Vibrates on Teeth/Hard Surface



32-40KHZ

Vibration frequency



1500mAH

Long Battery Life



< 60db

Noise index



Start Vibrating
touch the teeth



Stop Vibrating
touch the soft gums

SAFE TO THE GUMS



Touch with hard teeth

✓ Start working



Touch with soft gums

✗ Stop working



LED Light to illuminate oral situation

10 Min Timer to Auto-Stop

NEGATIVE ION GENERATOR EFFICIENCY IMPROVEMENT BY 50%

Grasp the ion connection ring to generate ions
Negative ions are transported to teeth and
repelled by negative ions The principle of
removing dirt from teeth



ION CONNECTING RING



Chengdu Chengbei Technology Co., Ltd.



+852 9506 3884



candice@pet-belief.com



product-pet.com

Room 502A, Building 1, 38 Tianyi Street, High-Tech Zone, Chengdu, Sichuan, China